


SOLAWAVE

ADVANCED SKINCARE WAND
with Red Light Therapy

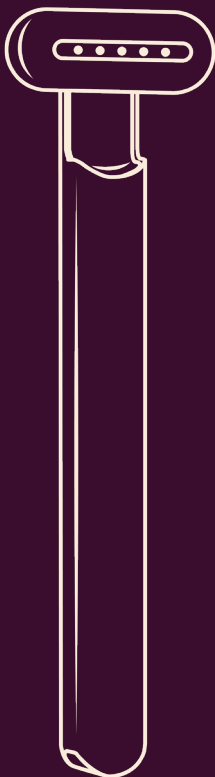

Scan to register your product
warranty and watch tutorials

ENGLISH INSTRUCTION MANUAL

ITEMS INCLUDED



Solawave Wand




Charging Cable



Instruction Manual

PRODUCT DIAGRAM

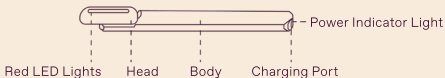

FEATURES



Microcurrent

Microcurrent helps promote a refreshed appearance of the skin and neck.



Red Light Therapy

The Solawave Wand uses visible red light with a centroid wavelength of 660 nanometers. Red light therapy can deeply rejuvenate the skin.



Vibrating Massage

Facial massage can reduce the appearance of puffiness and tired skin.



Therapeutic Warmth

Warming the skin can temporarily relieve eye fatigue and lead to a beautiful, healthy glow.



Smart Touch Sensor

The Solawave Wand activates when it is held firmly against hydrated skin. The wand automatically deactivates when it is pulled away from the skin.



90° Rotating Head

Solawave's 90° rotating head allows the device to better adapt to the contours of your face. Rotate your Wand 90° clockwise to the horizontal position. To return the Wand to the vertical position, turn the wand 90° counterclockwise.

GETTING STARTED

- 1) Peel off the protective stickers on the base and head of the Solawave Wand.
- 2) Charge the Solawave Wand fully prior to its first use.

CHARGING INSTRUCTIONS

To charge, plug the smaller end of the provided charging cable into the USB port on the bottom of the Solawave Wand, and the larger end of the cable into a USB power source (such as a wall charger, car charger, computer or portable battery pack). Ensure the power source output is 5V/1A (or 5 watts). Do not use a 'fast charger' with the Solawave Wand. The indicator light at the bottom of the Solawave Wand will flash slowly to indicate that charging is in progress. The product takes up to two hours to fully charge.

POWER INDICATOR

- When the battery is low, the indicator light will flicker.
- While charging, the indicator light will blink on and off.
- When fully charged and still connected to a power source, the indicator light will stop blinking and remain on.
- When fully charged, the battery will provide up to 60 minutes of usage.

HOW TO USE

Tip: visit [solawavetutorials.com](https://www.solawave.com/tutorials) to watch instructional tutorials that demonstrate how to use the Solawave Wand.

- 1) Clean your face and remove makeup. Apply a Solawave Activating Serum to your skin (or a serum or moisturizer of your choice). The Wand will not reliably activate if your skin is not properly hydrated.
- 2) Apply the Solawave Wand to the desired treatment areas designated in the 'Getting Started' manual included with your Wand.
- 3) We recommend rotating the head of the Solawave Wand 90° when gliding it over your face and neck area. The eye area is best treated with the Wand in its original vertical position. To return the Wand to its original vertical position, rotate the Wand 90°. See the diagram below to ensure proper usage of the Wand:



RECOMMENDED USAGE

We recommend starting out at 5 minutes per day for the first few days. As your skin acclimates, you can use it for longer and as frequently as desired for up to 30 minutes per day. Always be mindful of how your skin reacts to beauty devices and decrease usage if redness or irritation occurs.

IMPORTANT

- In order to preserve battery health, we recommend unplugging the Solawave Wand once it's finished charging.
- Do not disassemble or modify this product.
- Follow proper instructions to rotate the head. Do not rotate the head in the incorrect direction or it will break
- Do not rinse directly with water. Instead, wipe with a lightly dampened cloth after each use. For additional sterilization, spray the cloth with an alcohol-based sanitizer prior to wiping the device.
- Do not leave the Wand out in direct sunlight for extended periods of time.
- Do not activate or use the Wand while it is connected to a power source, this will decrease battery life.
- Keep this product out of the reach of children.
- Save your purchase receipt in case you need to file a warranty claim.

PRODUCT PROTECTION PLAN

We recommend purchasing our Protection Plan, which provides a complementary replacement in the event of accidental damage (e.g. broken head). Customers can purchase this plan up to one month from the date of purchase. To sign up and view full terms, visit solawaveplan.com

LIMITED WARRANTY

Solawave warrants that this product is free of defects in materials and workmanship for one (1) year from the date of purchase. You must register your Solawave purchase in order to be eligible for the limited warranty. Visit solawavewarranty.com to register your purchase. This warranty is valid only in accordance with the conditions set forth below:

- Normal wear and tear is not covered by this warranty
- This one (1) year warranty applies to consumer use only and is void when the product is used in a commercial or institutional setting.
- This warranty extends only to the original consumer purchaser and is not transferable. In addition, proof of purchase must be demonstrated. This warranty is void if the product has been subject to accident, misuse,

abuse, improper maintenance or repair, unauthorized modification, or any other use not found in the printed instructions.

This limited warranty is the only written or expressed warranty given by Solawave. Any implied warranties on the product (including but not limited to any implied warranties of merchantability or fitness for a particular purpose) are limited in duration to the duration of this warranty. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Repair or replacement of the product is the exclusive remedy of the consumer under this warranty. Solawave shall not be liable for any incidental or consequential damages for breach of this warranty or any implied warranty on this product. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

PROCEDURE FOR WARRANTY REPAIRS OR REPLACEMENT:

If warranty service is necessary, the original purchaser must contact Solawave for warranty service instructions at support@mysolawave.com.

For your protection, it is suggested you insure the shipment of the device back to Solawave. Damage that occurs during shipment is not covered. Packages that are missing proper documentation and accessories will not be processed.

INDICATIONS FOR USE

Over-the-counter cosmetic use on the face and neck.

CONTRAINDICATIONS FOR USE

The Solawave Wand should not be used in the following areas:

- Breast area
- Larynx
- Chest
- Groin area

WARNINGS

Do not submerge the unit or charging cable in water or any other liquid. If you store your Solawave Wand and/or charging cable in an area such as a bathroom counter, ensure the product does not come in contact with water or other liquid. If the base of the Wand or its head is exposed to water, it will stop working properly.

THIS SOLAWAVE WAND SHOULD NOT BE USED BY:

- People under the age of 18
- Pregnant women
- People subject to seizures
- People with skin cancer
- People with a cardiac pacemaker
- People with implanted defibrillators/stimulators
- People with implanted electronic devices

If you have any medical concerns, medical illness, medical conditions such as epilepsy or seizures, or recent facial surgery please consult your doctor before using the Solawave Wand.

Electronic monitoring equipment (such as ECG monitors and ECG alarms) may not operate properly when the Solawave Wand is in use. The Solawave Wand should only be used over healthy, intact skin.

The Wand should not be over, or in proximity to cancerous lesions, or applied over swollen, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc).

The Solawave Wand is not designed to activate when in contact with facial hair. We recommend shaving before use as hair can interfere with the conductivity of the wand. The facial area under a beard or mustache cannot be treated.

Potential adverse reactions may include skin irritation or redness. For adult use only. Keep out of reach of children.

Do not use the wand across your chest as the introduction of electrical current into the chest may cause rhythm disturbances to your heart.

PRECAUTIONS

The Solawave Wand is designed for cosmetic use only and for individuals in good health.

A slight tingling sensation may occur while using the Solawave Wand. There also may be a perception of flashing lights during the use of the wand. If you notice this condition constantly and/or without the use of your Solawave Wand, consult your doctor.

Please review the user manual and online instructional videos at solawavetutorials.com before using your Solawave Wand.

Caution should be exercised for patients with suspected or diagnosed heart problems.

Caution should be exercised for patients with suspected or diagnosed epilepsy.

Caution should be exercised if you have a tendency to hemorrhage following acute trauma, and/or when using the Solawave Wand over areas of the skin that lack normal sensation.

If you experience skin irritation or hypersensitivity due to the electrical stimulation or serum (e.g., the Solawave Activating Serum), you can likely reduce the irritation by using an alternate serum/gel.

The Solawave Wand should not be used while driving, operating machinery, or during any activity in which involuntary muscle contractions may put the user at undue risk of injury.

If you have suspected or diagnosed heart disease, you should follow precautions recommended by your physician. Consult your physician prior to using the device after a recent surgical procedure.

ADVERSE REACTIONS

If you experience a persistent headache, skin irritation, or other adverse reactions during or following the use of the Wand, you should stop using the device and consult your physician.

VERY IMPORTANT

This device is intended for cosmetic use on the face and neck. Any harmful consequence resulting from misuse, or application to other body areas, connection to improper voltage sources, application of 3rd party serums, or any other improper applications is neither the responsibility of Solawave nor its affiliates.

TAG US

We'd love to see your progress as you incorporate Solawave into your beauty routine! Share your photos and videos with the community by tagging @solawave on Instagram™.

PRODUCT SPECIFICATIONS

NAME: Advanced Skincare Wand
with Red Light Therapy

MODEL NO.: SW001, SW002, SW003, SW012

MATERIAL: Aluminum Alloy & PC

SIZE: 0.5" W x 0.6" H x 5.8" L

NET WEIGHT: 1 oz.

BATTERY CAPACITY: 350mAh

CHARGING TIME: 2 hours

USAGE TIME: 1 hour

INPUT VOLTAGE: 5V

POWER: 1.5W

CERTIFICATIONS: RoHS



CONTACT

support@mysolawave.com

For more information, visit [mysolawave.com](https://www.mysolawave.com)

DISTRIBUTED BY

UVVU Inc
2633 Lincoln Blvd. #1025
Santa Monica, CA 90291



UVVU Inc
2633 Lincoln Blvd. #1025
Santa Monica, CA 90291

EU AUTHORIZED REPRESENTATIVE

Obelis s.a.

Bd. Général Wahis 53 B-1030, Brussels, Belgium